

# SpectraMax i3x Multi-Mode Detection Platform

Explore a wealth of applications in one future-ready system



## Benefits

- User-upgradeable application modules including cellular imaging
- Sensitivity across spectrum with Spectral Fusion™ Illumination
- Expanded dynamic range with cooled PMT
- Control and analytics provided by SoftMax® Pro Software



# SpectraMax i3x


The SpectraMax i3x from Molecular Devices is a multi-mode detection system that evolves with your future needs and offers an unlimited breadth of application possibilities.

## Future ready

The SpectraMax® i3x Multi-Mode microplate reader measures spectral-based Absorbance, Fluorescence, and Luminescence with the added functionality of modular upgrades for Western Blot, Imaging, and Fast Kinetics with Injectors.

The SpectraMax i3x reader allows you to unravel the mysteries of science by exploring cellular pathways and protein activation and expression in one system. Protect your initial investment — purchase a system with the flexibility to add novel detection capabilities without the need for service engineers or costly system downtime. The SpectraMax i3x reader grows with you as your research areas expand.

Be Future Ready with the SpectraMax i3x Multi-Mode Detection Platform.

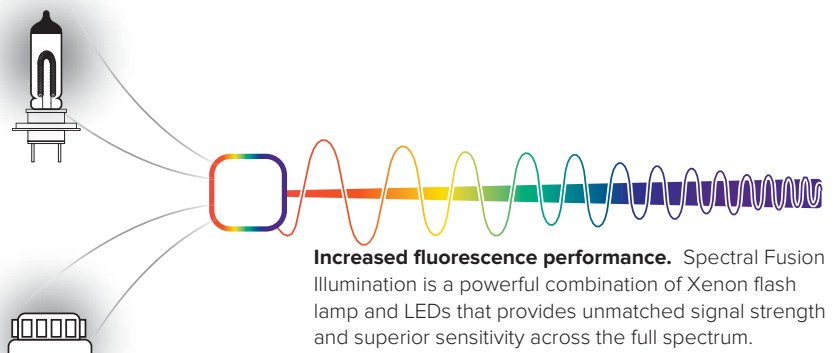

**Increased fluorescence performance.** Spectral Fusion Illumination is a powerful combination of Xenon flash lamp and LEDs that provides unmatched signal strength and superior sensitivity across the full spectrum.

## ENGINEERED TO PERFORM

### Expanded dynamic range

Engineered for performance with spectral fusion illumination for increased sensitivity across the entire excitation range and a cooled photomultiplier tube (PMT) for improved detection in extremely low light. Generate more data points without the need to dilute.

### User-upgradeable applications


User-exchangeable detection modules expand the system's detection capabilities to include Time-Resolved Fluorescence, HTRF, Fast Kinetics with Injectors, and Western Blot detection.

Future proof your discoveries and let the SpectraMax i3x reader evolve with your research.

### One complete solution

With available options such as the SpectraMax® MiniMax™ 300 Imaging Cytometer, ScanLater™ Western Blot cartridge, reagents optimized for high performance, and the industry-leading data acquisition and analysis tool SoftMax Pro, the SpectraMax i3x Detection Platform is the total solution for all your research needs.

Let Molecular Devices® and the SpectraMax i3x Multi-Mode Detection Platform power your next landmark discovery.

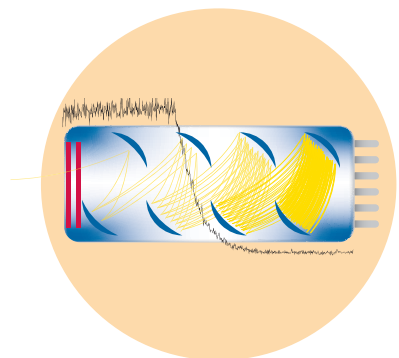

### Quantitative low-light measurement

Cooled PMT reduces background noise allowing for a more sensitive, wide dynamic range in extremely low light.



### New applications in minutes

Adding modes and functionality is just minutes away. Simply insert a cartridge to expand your application capabilities.



### Your total solution

Matched reagents, an imaging cytometer option, user-upgradeable detection modules, and industry-leading SoftMax Pro provide you with the ability to explore every pathway.



## A WEALTH OF APPLICATIONS IN ONE FUTURE-READY SYSTEM


### Investigate every aspect of a cellular pathway

From imaging of cell confluence and viability under different treatment conditions to quantitation of nucleic acids and protein to western blot analysis, a wealth of new knowledge is captured using a single instrument. One versatile software package powers data acquisition and analysis, from raw data to publishable results. See more and do more cellular to protein analysis using a single detection system.




### Live cell imaging

The SpectraMax MiniMax 300 Imaging Cytometer option allows for live cell images and analysis (below, left). The ScanLater Western Blot Detection Cartridge enables protein detection (bottom, right).

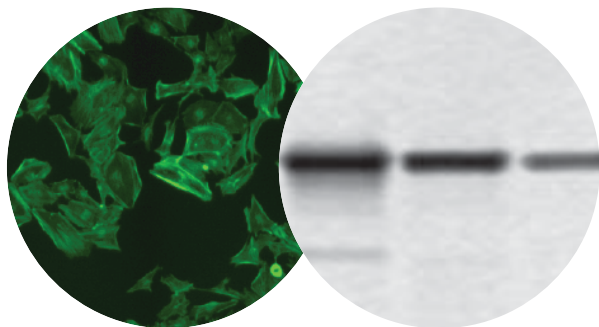

### Visualize cells with your microplate reader



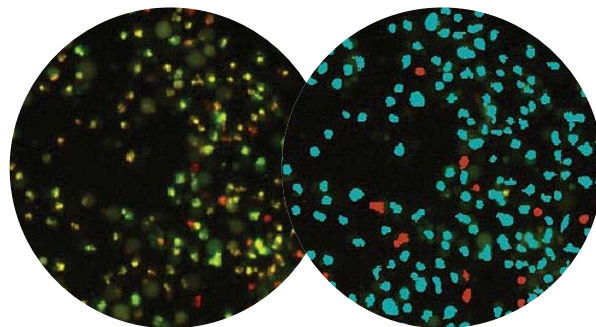




Imaging with the SpectraMax MiniMax 300 Imaging Cytometer mirrors the plate reading workflow on the SpectraMax i3x System. The plate is set up for reading and images are acquired according to specified parameters. Cells in each image are identified by SoftMax Pro Software and cell-by-cell statistics are collected. Data are then analyzed and visualized in different graphical representations.




### Optional enhancements

- SpectraMax MiniMax 300 Imaging Cytometer
- ScanLater Western Blot System
- SpectraDrop™ Micro-Volume Microplate
- SpectraTest® Validation Packages (ABS1, FL1, LM1)
- SpectraMax i3x Injector Cartridge
- SoftMax Pro GxP Microplate Data Compliance Software
- IQ/OQ Protocols
- Additional detection modules

# Additional cartridges

Cartridge	Description	Name	Specifications	Optimized sensitivity	Guaranteed sensitivity	Slots used
 ScanLater	Western blot detection using ScanLater Western Blot Assay Kit TRF-based with 340/80 nm EX and 616/10 nm EM	0200-7027	EX range: 340/80 nm EM range: 616/10 nm	High fg levels of Streptavidin	High fg levels of Streptavidin	2
 AlphaScreen	AlphaScreen and AlphaLisa detection using 1 W 680 nm EX laser diode and a 570 nm (100) EM filter	0200-7017POS	Alpha 384 STD 96- and 384-well plates		< 100 amol (384-well)	1
	Pick best speed, sensitivity, and price for your needs	0200-7018POS	Alpha 384 HTS 96- and 384-well plates		< 100 amol (384-well)	1
	Guaranteed sensitivity: < 100 amol phosphorylates biotin-peptide in 25 µL assay volume in a 384-well plate	0200-7019POS	Alpha 1536 HTS 96, 384, and 1536-well plates		< 100 amol (384-well)	1
 HTRF	Cisbio HTRF detection with optimized Xenon light source and 616, 665 nm EM filters Measures both emissions simultaneously	0200-7011POS	6- to 1536-well plates	Exceeds Cisbio certification requirements		2
 TRF	LED light source and Europium EX and EM filters (370-616 nm) Suitable for assays using Europium chelate and similar labels Includes 642 nm EM filter for TR-FRET assays with Samarium labels	0200-7008POS	6- to 1536-well plates	96–0.03 pM 384–0.03 pM 1536–0.125 pM	96-well: 0.1 pM 384-well: 0.1 pM 1536-well: 0.375 pM	1
 FP	Fluorescence Polarization detection for fluorescein- or rhodamine-like labels Using specific LED and EX/EM filters for 6- to 1536-well plates	0200-7009POS 0200-7010POS	Fluorescein FP EX 485 nm, EM 535P and 535S nm Rhodamine FP EX 535 nm, EM 595P and 595S nm	96–1.0 mP 384–1.5 mP 1536–2.0 mP	96-well: 3 mP 384-well: 4 mP 1536-well: 6 mP	1
 Injectors	Dual Auto Injectors to expand your research capabilities to include flash reactions	Available October 1, 2015				
 Custom Solutions	Custom cartridges are available and designed to meet your specific research needs	Ask your local sales representative for more information				

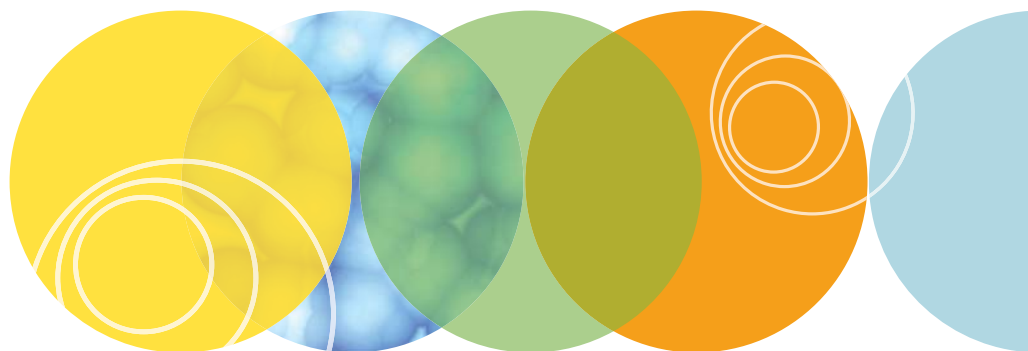



# Technical specifications (base system)

General specifications	
Dimensions (in.)	12.63 (H) x 15.38 (W) x 23.38 (D)
Dimensions (cm)	42.23 (H) x 39.05 (W) x 59.37 (D)
Weight	68.3 lbs. (31.0 kg)
Power consumption	< 200 watts
Power source	100–240 Vac, 2 A, 50/60 Hz
Robotic-compatible	Yes

General photometric performance	
Plate formats	6 to 1536 wells <sup>§</sup>
Light source	Spectral Fusion Illumination (Xenon flash lamp + high-powered LEDs or laser diode in detection cartridges)
Reading capabilities	Microplates, cuvettes (via adapter)
Detectors	PMT and/or photodiode
Shaking	Linear and orbital
Dual Auto Injectors	Available October 1, 2015
Temp. control	4°C above ambient to 45°C
Temp. uniformity	± 0.75°C
Temp. accuracy	±1°C at 37°C set point
Environmental control	Gas quick connect
Spectral scanning	Abs, FI, Lum
Endpoint reading	All modes
Kinetic reading	All modes
Well scanning	Over 20 by 20 in all modes
Wavelength selection	1.0 nm increments

Standard read times (minutes:seconds)*		
	96 wells	384 wells
Absorbance	0:30	1:40
Fluorescence intensity	0:25	1:25
Luminescence	0:30	1:15

Absorbance photometric performance	
Wavelength range	230–1000 nm
Wavelength bandwidth	4.0 nm
Wavelength accuracy	± 2.0 nm
Wavelength repeatability	± 1.0 nm
Photometric range	0–4.0 OD
Photometric resolution	0.001 OD
Photometric accuracy	< ±0.010 OD ±1.0%, 0–2 OD
Photometric precision	< ±0.003 OD ±1.0%, 0–2 OD
Stray light	< 0.05% @ 230 nm

Fluorescence intensity performance		
Wavelength range	250–850 nm	
Wavelength selection	1.0 nm increments	
Bandwidth (EX/EM)	Adjustable EX 9/15 nm EM15/25 nm	
Dynamic range	> 6 logs	
<b>Top sensitivity (fluorescein)</b>	<b>Optimized</b>	<b>Guaranteed**</b>
96 wells	0.5 pM	3 pM
384 wells	1 pM	4 pM
<b>Bottom sensitivity (fluorescein)</b>	<b>Optimized</b>	<b>Guaranteed**</b>
96 wells	5 pM	10 pM
384 wells	5 pM	20 pM

Luminescence performance		
Wavelength range	300–850 nm	
Wavelength selection	Choice of simultaneous detection of all wavelengths or selection in 1.0 nm increments	
Dynamic range	> 6 decades	
Cross-talk	< 0.4% in white 96- and 384-well microplates	
<b>Sensitivity (ATP-Glow)</b>	<b>Optimized</b>	<b>Guaranteed**</b>
96 wells	3 pM	15 pM
384 wells	6 pM	30 pM

Scan the QR code to request more information on the SpectraMax i3x Multi-Mode Detection Platform.



\* With 6 flashes in absorbance and 3 flashes in fluorescence mode and 0.1 sec./well integration in 96-well luminescence mode and 0.04 sec./well integration in 384-well luminescence mode

§ 1536 detection available via detection cartridges

\*\* For properly functioning, operating, and maintained equipment

## Contact Us

Phone: [800.635.5577](tel:800.635.5577)  
 Web: [www.moleculardevices.com](http://www.moleculardevices.com)  
 Email: [info@moldev.com](mailto:info@moldev.com)

Check our website for a current listing of worldwide distributors

## Regional Offices

USA and Canada  
 +1.800.635.5577  
 Brazil  
 +55.11.3616.6607  
 China (Beijing)  
 +86.10.6410.8669

China (Shanghai)  
 +86.21.3372.1088  
 Germany  
 00800.665.32860  
 Japan (Osaka)  
 +81.6.7174.8831

Japan (Tokyo)  
 +81.3.6362.5260  
 South Korea  
 +82.2.3471.9531  
 United Kingdom  
 +44.118.944.8000